


# Free2Go Rollator

## Combining Bathroom Safety and Mobility

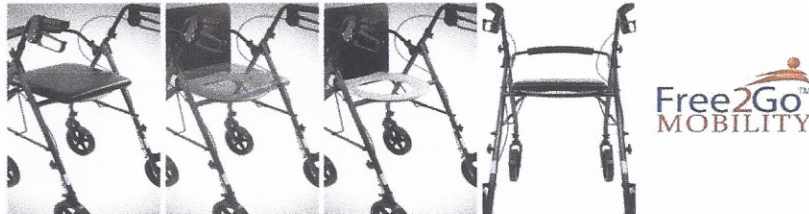

The Free2Go Rollator combines the benefits of the rolling walker with a toilet safety product, thereby allowing one to safely, confidently and discreetly use a toilet in *any* restroom - at home, in the homes of others, and in public restrooms.

Upon entering a restroom, simply push the rollator over a standard toilet, engage the wheel locks, and raise the padded seat to use. The flip-up padded seat is equipped with a pouch to conveniently store disposable wipes to clean the unit after each use.

- Combines the benefits of a raised toilet seat, toilet safety frame and rolling walker
- The frame of the rollator provides support when sitting, lifting and standing
- Provides additional height over a standard toilet to reduce body strain
- Restores dignity and independence
- Reduces the fear of falling and the possibility of life-threatening accidents
- Portable, lightweight and discreet design
- 8" flat wheels for added stability
- Removable flip-up curved padded backrest
- Vinyl padded seat with easy-folding system
- Storage pouch located under seat to store cleaning wipes
- Limited lifetime warranty on frame and 90 days on parts

### Product Specifications:

Seat Measurements: 15" (W) x 15 ½" (D)

Toilet Seat Opening: 8 ¼" x 9 ½"

Seat-to-Floor Height: 23 ½" (top of seat)

Seat-to-Floor Height: 22 ¼" (top of toilet seat)

Adjustable Handle Height: 33 ½" – 38"

Unit Width: 24"

Product Weight: 20.2 lbs.

Patient Weight Capacity: 300 lbs.